

# CURADLE' PET BROODER



CURADLE



세계일류상품

**CURADLE offers the best  
in luxury and practicality for your business.**

Available for optimal environment settings depending on pet's condition

# Nursery Hospital Brooder

: The highest level of Pet Brooder

## ANION ...

The reason why we feel refreshing when we are in woods or mountain is that there are much anion, the atom element with electric charge. If we get anion, it increase the vitality by accelerating metabolism of cell. And, it cleans our blood and helps recovering from fatigue, heavy sleep, and improvement of appetite.



## MAIN FUNCTIONS

- ▷ Innovatively improve stability and durability with centralized heating system applying PTC memory element from carbon filmed heating system
- ▷ Built-in air filter which removes indoor dust
- ▷ Built in Anion activation function
- ▷ Water pump included for humidification
- ▷ Reduce up to 40% of electric power consumption compare with the existing products by using PTC heater

## EASY FUNCTIONS

- ▷ Automatic Temp/Humidity control and setting
- ▷ Recessed brooding tray for easy cleaning
- ▷ 5 stage dimmer device for indoor illumination control
- ▷ Rotary Door
- ▷ PTC heated humidification system kills germ
- ▷ Reduce of fan noise inside of brooders and provide better pleasant environment (69dB)


## OTHER FUNCTIONS

- ▷ Nubulizer(Respiratory therapy apparatus) and Oxygen tank(Oxygen therapy apparatus)
- ▷ BLDC fan for optimum control of inner conditions
- ▷ Elegant outer design and optimum inner space distribution
- ▷ Alarm function in case of abnormal temperature
- ▷ It is possible to stack machine which enable to provide space efficiency (Brooder standing - Option)





OPERATIONAL FUNCTIONS

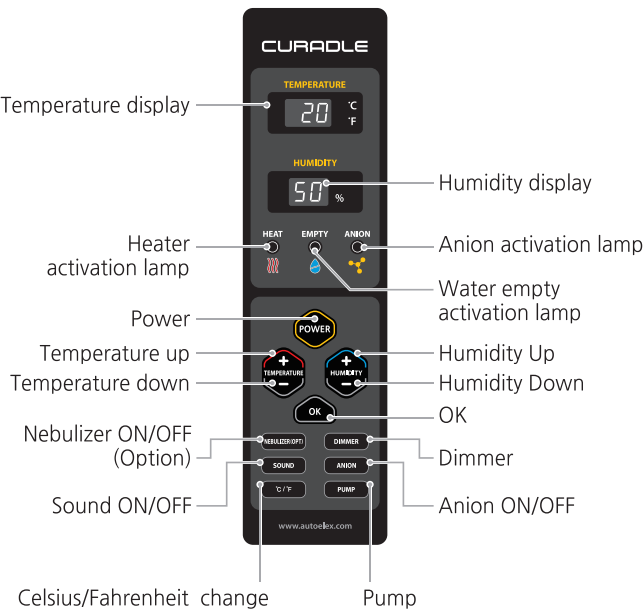

- You can control the optimum conditions according to the state of pet health

Classify	Normal Use	In case of sickness
Temperature	OFF ~ 32[°C] 89[°F]	25 ~ 36[°C] 77 ~ 96[°F]
Humidity	OFF ~ 50[%]	50 ~ 60[%]
Anion	ON or OFF	ON or OFF
Other options		Nebulizer / Oxygen Tank

※ It can be different according to the species of pets.

Connecting Oxygen & Auxilioary equipment

- Nubulizer and Oxygen tank can be connected and used separately or simultaneously.



This is Nebulizer type. Nebulizer device & Oxygen tank sold separately.

Installation of secondary tray

- Secondary tray can be installed and used as a litter tray for cats or a birth box, etc.



※ The secondary tray is sold separately.


INTERIOR



Water cap &  
Water nipple



Nebulizer & Oxygen  
connecting part

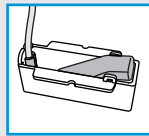

Humidification  
unit



PTC Heater



BLDC Fan



Water pump




LED Indoor light



Micro processor  
controller



Brooder standing  
(Sold separately)

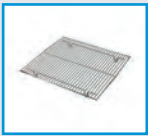

Stainless mesh  
(Sold separately)




Secondary tray  
(Sold separately)




Nebulizer  
(Sold separately)

● PET BROODER has two different sizes.

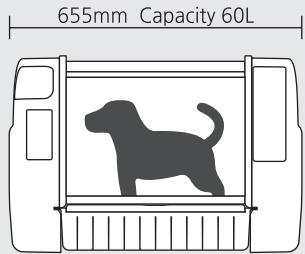

[PET BROODER 60]

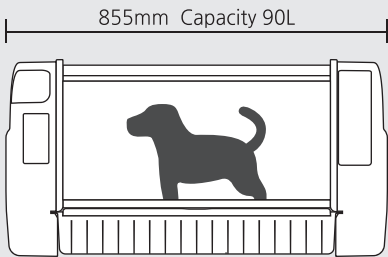

[PET BROODER 90]

SPECIFICATION

Product name	CURADLE PET BROODER 60	CURADLE PET BROODER 90
Model No.	MX-B60N	MX-B90N
Power supply voltage	AC 100~120[V], 50/60[Hz] / AC 220~240[V], 50/60[Hz]	
Power consumption	Average: 80[W]	Average: 90[W]
Temperature control range	Ambient temperature ~ 38[°C]   Ambient temperature ~ 100[°F]	
External environment temperature	Optimum range : 20 ~ 25[°C]   68 ~ 77[°F]	
Humidity control range	40 ~ 60[%] (It can be changed according to the conditions)	
External environmental humidity	40 ~ 60[%]	
Weight	N.W 9[Kg]   G.W 13[Kg]	N.W 10[Kg]   G.W 14.5[Kg]
Size	(W)655 x (L)470 x (H)440[mm]	(W)855 x (L)470 x (H)440[mm]
Fuse Standard	250[V] 10.0[A] (Ø5*20[mm])	

※ To improve the performance of the product, specifications are subject to change without prior notice.